

# GlassOuse

Assistive Device


[www.glassouse.com](http://www.glassouse.com)

**GS05 Press Switch 100g**

User Manual

## Press Switch 100g GS05

Dimensions	Diameter: 96.9mm x Thickness: 30.55mm
Weight	68±2g
Cable Length	185±5cm
Activation Surface	86.30mm
Activation Force	100g

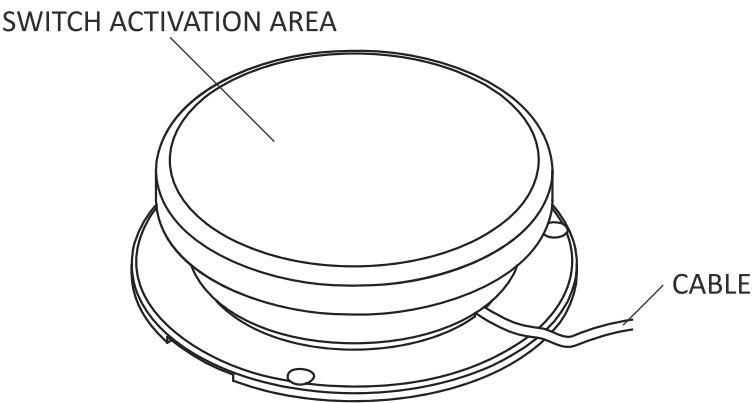

## Press Switch 100g

A durable switch option for users with less physical control. GS05 Press Switch provides 8.63cm activation surface, making the item comfortable enough for both user groups needing either a large or smaller target area to interact. Operated by 100g activation force, GS05 Press Switch has an auditory click and tactile feedback.

Press Switch comes with a 185cm long cable and a 3.5mm mono jack plug-in. Can be used with different mounting options or double-sided tapes.

## SWITCH CONNECTION

---

1. Plug GS05 Press Switch to GlassOuse 3.5mm jack adapted device via mono 3.5mm jack connector. 2. Plug GS05 press switch to Glassouse 3.5mm Jack adapted device via Mono 3.5mm Jack Connector.
  - 2.1 Place GS05 Press Switch on a flat surface and adjust the position until you feel comfortable enough to do the pressing action before using double-sided tapes to secure the switch.
- OR
- 2.2 Use a mounting option to secure the position of GS05 Press Switch.

## WARNINGS

---

While removing your Press Switch from your GlassOuse or other devices, pay attention not to damage the wire joining points on the 3.5mm jack connector and do not pull the cable harshly.