

GlassOuse

Assistive Device


www.glassouse.com

GS08 Touch Switch

User Manual

TOUCH SWITCH GS08

Dimensions	Diameter: 96.9mm x Thickness: 28.55mm
Weight	125±2g
Cable Length	185±5cm
Battery Type	3x AAA Batteries

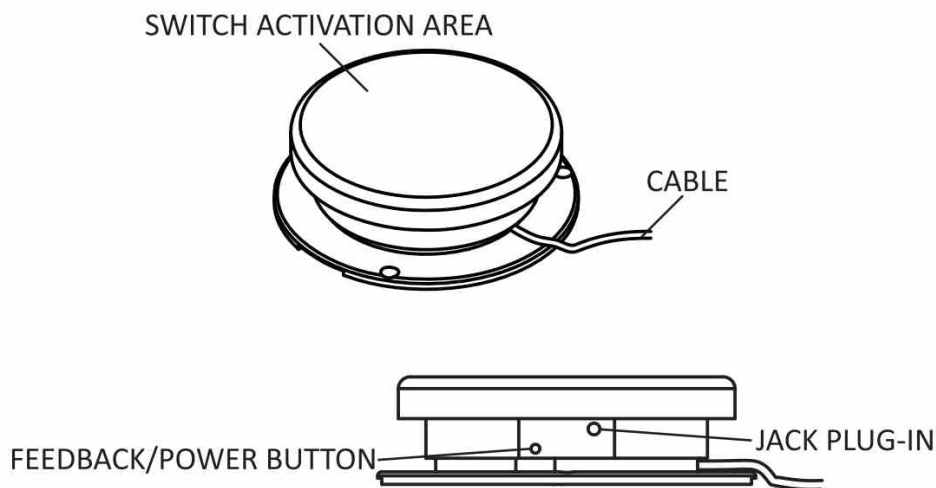

Touch Switch

With an 8.63cm activation surface, Touch Switch is highly sensitive to detect slight touches on activation area without any delay. Simply touch any area on the switch module to perform your clicks including the edge and sides of the device GS08 Touch Switch comes with a 3.5mm mono jack connector plug and it is compatible with any kind of adapted devices besides GlassOuse.

Touch Switch works with 3x AAA batteries and its smart battery system turns the device off in one hour if it is not used. Auditory and visual feedback can be selected according to the user's needs easily or can be turned off if desired. With the help of the external secondary switch, you can extend the using range of Touch Switch with a smaller activation module.

FEEDBACK SIGNALS

You can set the click feedback according to your needs. Use the tiny power button on the switch module for the desired feedback type as following. A second of press is enough to change the feedback mode.

1. Press the power button once - only visual LED feedback on without any auditory sound.
2. Press the power button once again - only auditory sound feedback on without no LED.
3. Press the power button once again - both visual LED and auditory sound feedback on together.
4. Press the power button once again - no visual LED or auditory sound feedback.

SWITCH CONNECTION

1. Plug Touch Switch to GlassOuse or any 3.5mm jack adapted device via 3.5mm connector.
2. Adjust the position of Touch Switch until you feel comfortable enough to touch the switch.
3. Press the power button on your switch for 3 seconds to activate Touch Switch. The LED light will flash with beep sound. Press the power button again for 3 seconds to turnoff your Touch Switch. The LED light will flash with beep to indicate that the switch is turned off.

WARNINGS

1. While removing your Touch Switch from your GlassOuse or other devices, pay attention not to damage the wire joining points both on the 3.5mm jack connector and the Touch Switch module parts.
2. To change the battery, remove the tiny screws at the bottom. Pay attention to the + and - pole remarks of the battery box during installation.
3. Keep Touch Switch away from magnets.

EXTERNAL SECONDARY SWITCH

You can combine your Touch Switch with a tiny and more portable option to extend the using range. Simply plug the secondary switch and place it according to your needs.